

Primerdesign

R01050

Norovirus Genogroups I and II

Kit version: 2

Target region:
RNA polymerase gene

genesig[®] Easy RNA Kit

for use on the genesig[®] q16
50 reactions

GENESIG

Kits by Primerdesign

For general laboratory and research use only

Norovirus Genogroups I and II

Norovirus is known to cause acute gastroenteritis. It is a small (27-38 nm), round, non-enveloped RNA virus belonging to the Caliciviridae family and is responsible for over 80% of non-bacterial outbreaks of gastroenteritis in the world. It affects individuals of all ages, with a distinct seasonal link to winter. It has a genome of 7.6 kb that is positive sense and has a single stranded linear confirmation. It encodes a major structural protein (VP1) of about 58 to 60 kDa and a minor capsid protein (VP2).

Transmission occurs predominantly through ingestion of contaminated water, food and airborne transmission, as well as contact with contaminated surfaces. The ease with which norovirus is transmitted and the low infectious dose required to establish an infection results in extensive outbreaks in numerous environments, such as hospitals, hotels and schools. There is no antiviral drug available to treat this infection and little is known about its pathogenicity. However, it has been observed that the virus can be taken up by enterocytes where translation of viral non-structural proteins can occur; it damages and alters intestinal microvilli, leaving them blunt and broadened, thus inhibiting absorption; it causes crypt cell hyperplasia and leads to apoptosis of enterocytes.

An incubation period of 24-48 hours is usual. Infection is characterized by the acute onset of nausea, vomiting, abdominal cramps, aching limbs, raised temperature and diarrhoea that generally last for about 48 hours. However, more severe and prolonged infection may be observed in children and the elderly. There are now ten recognised norovirus genogroups, of which three (GI, GII, and GIV) are known to affect humans. Of these three genogroups, GI and GII are responsible for the majority of human infections and, since 2002, variants of the GII.4 genotype have been the most common cause of norovirus outbreaks. There have been 51 different genotypes identified within the genogroups, with a wide degree of genetic variability present even within each genotype.

Specificity


This Norovirus Genogroups I and II kit is designed for the in vitro quantification of Norovirus GI and GII genomes. The kit is designed to have a broad detection profile. Specifically, the Norovirus GI and GII assays will detect over 93% and 95%, respectively, of sequences available on the NCBI database at the time of last review.

The Norovirus GII assay is predicted to cross react with Noroviruses GIX, GVI and GIV. Newly emerging evidence suggests that there will be some cross reactivity with bacteriophages found in the human gut, however the effect of this is predicted to be negligible in samples from patients suffering from gastrointestinal infection.


The dynamics of genetic variation means that new sequence information may become available after the initial design. If you require further information or have a specific question about the detection profile of this kit, please send an email to techsupport@primerdesign.co.uk and our team will answer your question.

genesig[®] Easy: at a glance guide


For each RNA test

Component	Volume	Lab-in-a-box pipette	
Norovirus G1_v3.0 or Norovirus G2_v2.0 primer/probe mix	5 µl	●	
Your RNA sample	5 µl	●	
oasig [®] PLUS OneStep Master Mix	10 µl	●	











For each positive control

Component	Volume	Lab-in-a-box pipette	
Norovirus G1_v3.0 or Norovirus G2_v2.0 primer/probe mix	5 µl	●	
<u>Positive control template</u>	5 µl	●	
oasig [®] PLUS OneStep Master Mix	10 µl	●	

For each negative control

Component	Volume	Lab-in-a-box pipette	
Norovirus G1_v3.0 or Norovirus G2_v2.0 primer/probe mix	5 µl	●	
<u>RNase/DNase free water</u>	5 µl	●	
oasig [®] PLUS OneStep Master Mix	10 µl	●	

Kit Contents

Quantity	Component	Tube	Cap Colour
1	Norovirus G1_v3.0 primer/probe mix (50 reactions) FAM labelled for Norovirus G1; VIC labelled for the Internal Control		BROWN
1	Norovirus G2_v2.0 primer/probe mix (50 reactions) FAM labelled for Norovirus G2; VIC labelled for the Internal Control		YELLOW
1	Norovirus G1_v3.0 positive control template		RED
1	Norovirus G2_v2.0 positive control template		RED
1	Internal extraction control RNA		BLUE
2	oasig® PLUS OneStep Lyophilised 2X RT-qPCR Master Mix (50 reactions per glass vial)		GOLD
2	oasig® resuspension buffer For resuspension of the lyophilised Master Mix		ORANGE
1	RNase/DNase free water for resuspension of primer/probe mixes		WHITE
3	Template preparation buffer for resuspension of internal control template and positive control template		YELLOW
2	54 genesig® q16 reaction tubes		N/A

Reagents and equipment to be supplied by the user

genesig[®] q16 instrument

genesig[®] Easy Extraction Kit


This kit is designed to work well with all processes that yield high quality RNA and DNA but the genesig[®] Easy extraction method is recommended for ease of use.

genesig[®] Lab-In-A-Box

The genesig[®] Lab-In-A-Box contains all the pipettes, tips and racks that you will need to use a genesig[®] Easy kit. Alternatively, if you already have these components and equipment these can be used instead.


Step-by-step guide

1. Resuspend the test components

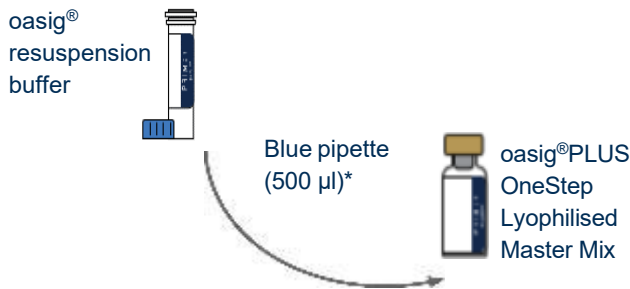

Use the blue pipette to transfer 500 µl* of the oasig® resuspension buffer into the tube of oasig®PLUS OneStep Lyophilised Master Mix and mix well by gently swirling.

*Transferring 525 µl of the oasig® resuspension buffer to your oasig®PLUS OneStep Lyophilised Master Mix (instead of the 500 µl recommended above) will enable you to take full advantage of the 50 reactions by accounting for volume losses during pipetting. In order to do so with the Lab-in-a-box fixed volume pipettes use 1x blue, 2x red and 1x grey pipettes to make the total volume. Please be assured that this will not adversely affect the efficiency of the test.

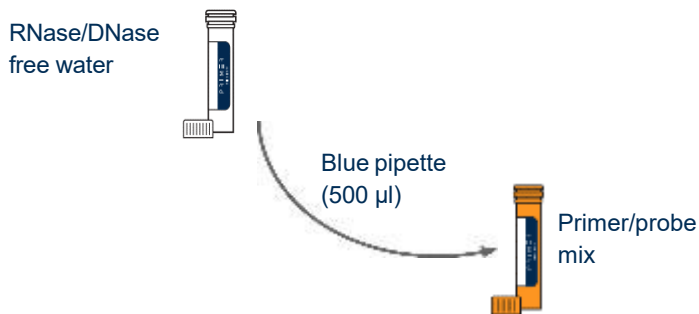


Then use the blue pipette to transfer 500 µl of water into the brown tube labelled oasig®PLUS OneStep Master Mix primer/probe mix. Cap and leave for 10 minutes to rehydrate. Then shake tube to mix. A thorough shake is essential to ensure that all components are resuspended. Use the same blue pipette to pipette up and down 10 times to allow complete resuspension. **Failure to mix well can produce poor kit performance.**

Store them in the freezer from hereon.

Top tip

- Ensure that the primer/probe mix is mixed thoroughly before each use by shaking and pipetting up and down 10 times.
- Once resuspended do not expose the genesig® Easy kit to temperatures above -20°C for longer than 30 minutes at a time.

2. Internal extraction control template

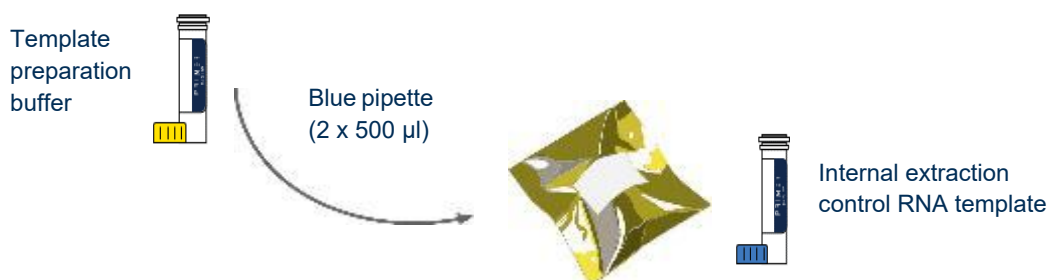

Use the blue pipette to transfer 1000 µl (2 x 500 µl) of template preparation buffer into the Internal Extraction Control RNA template tube. Cap and shake the tube to mix. Flick the tube to ensure that all liquid is collected at the bottom of the tube. Pipette up and down 10 times to allow complete resuspension.

The Internal Extraction Control RNA template should be added to your biological sample at the beginning of the RNA extraction process. It is extracted along with the RNA from your target of interest. The q16 will detect the presence of this Internal Extraction Control RNA template at the same time as your target. This is the ideal way to show that your RNA extraction process has been successful.


If you are using an extraction kit:

Use the red pipette to transfer 10 µl of Internal Extraction Control RNA template to your sample **after** the lysis buffer has been added, then follow the rest of the extraction protocol.

If using samples that have already been extracted:

Use the grey pipette to transfer 5 µl of Internal Extraction Control RNA template to your extracted sample.

3. Add primer/probe mix to all reaction tubes

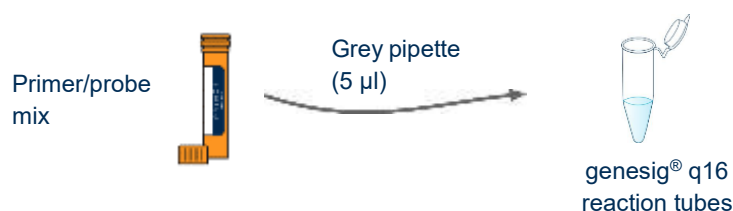


For every reaction to be run, use the grey pipette to add 5 µl of your Norovirus G1_v3.0 or Norovirus G2_v2.0 primers/probe mix to every genesig® q16 reaction tube.

Top tip

- Always pipette the primer/probe mix directly into the bottom of the tube.
- You can label the tube lids to aid your reaction setup but avoid labelling tube sides.

4. Add Master Mix to all reaction tubes

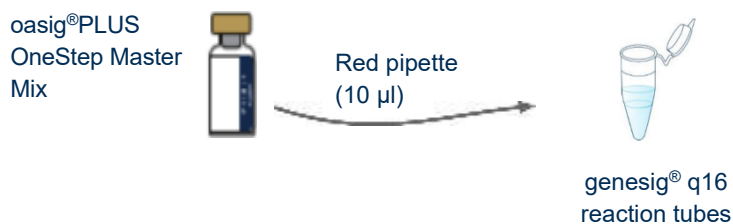


For every reaction to be run, use the red pipette to add 10 µl of the oasig®PLUS OneStep Master Mix to the reaction tubes containing primer/probe mix.

Move swiftly to begin your q16 run, as any delay after the oasig®PLUS OneStep Master Mix has been added can affect the sensitivity of your test.

Top tip

- Always add the oasig®PLUS OneStep Master Mix to the side of the tube to reduce the introduction of bubbles.

5. Negative control

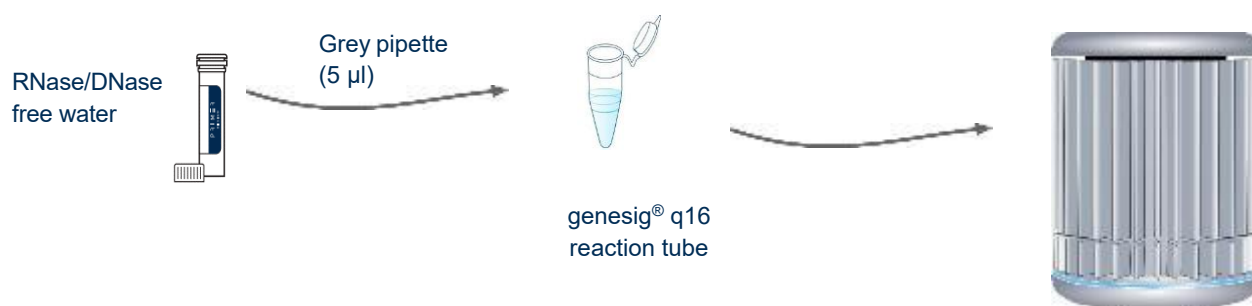

For each test you will require a negative control. Instead of RNA, water is used. The negative control sample should give a negative test result and thereby prove that any positive samples really are positive and not tested positive due to contamination.


To create a negative control reaction simply use the grey pipette to add 5 µl of the water to the required reaction tube (already containing primer/probe mix and oasig®PLUS OneStep Master Mix). Close this tube after adding the water.

Because some genesig® kit targets are common in the environment you may occasionally see a “late” signal in the negative control. The q16 software will take this into account accordingly.

Top tip

- Always add the water to the side of the tube to reduce the introduction of bubbles.

6. Set up a test

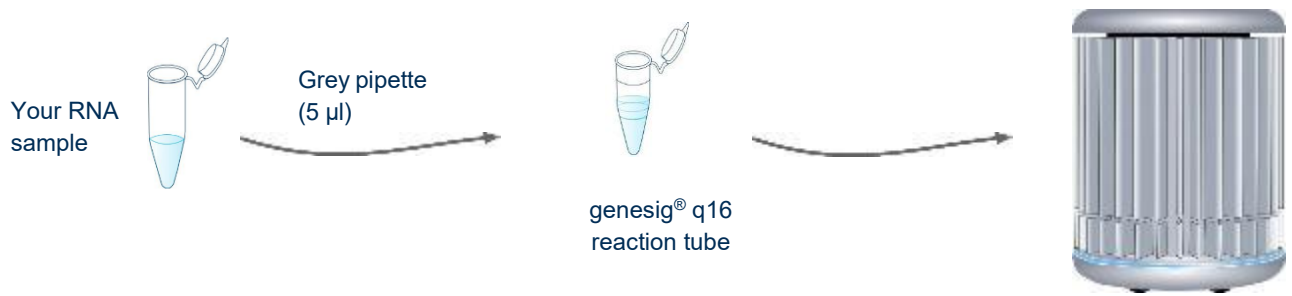


For each sample you wish to analyse, use the grey pipette to add 5 µl of your RNA sample to the required reaction tube (already containing primer/probe mix and oasis®PLUS OneStep Master Mix). Close these tubes after adding the sample. Always change pipette tips between samples.

Top tip

- Always add the RNA sample to the side of the tube to reduce the introduction of bubbles.

7. Positive control

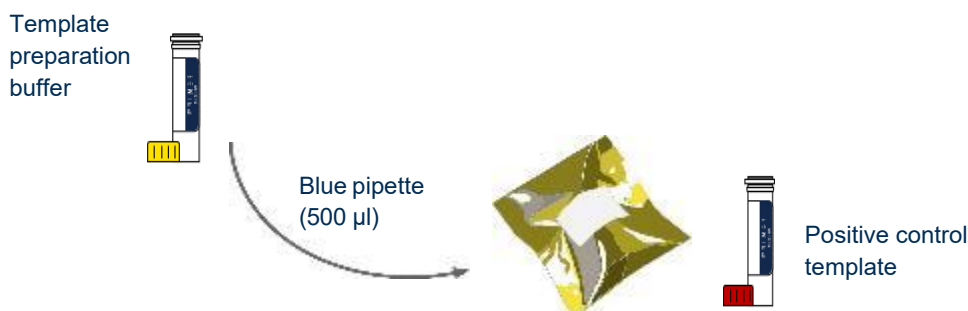


Use the blue pipette to transfer 500 µl of template preparation buffer into the positive control template tube. Cap and shake tube to mix.

Each time you run a test you will require a positive control. It serves two purposes:

1. It will always test positive so it shows that everything is working as it should be.
2. The q16 software knows how many copies of the target are present in the positive control. Therefore, it can automatically compare your sample of interest with the positive control to calculate the amount of target RNA in your sample.

To create a positive control reaction simply use 5 µl of the positive control template instead of your RNA sample and add this to the required reaction tube (already containing primer/probe mix and oasis® OneStep Master Mix). Close this tube after adding the positive control template.

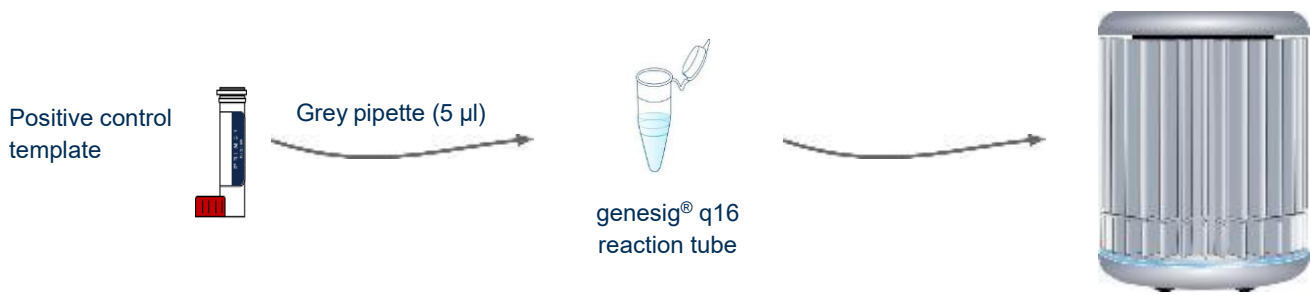


Take great care when setting up your positive control. The positive control template has the potential to give you a false positive signal in your other samples. Set positive controls up last after all other sample tubes are closed. Always change pipette tips between samples. You may even choose to set up positive controls in a separate room.

Top tip

- Always add the positive control template to the side of the tube to reduce the introduction of bubbles.

8. Running the test

Select the genesig Easy Target Detection Kit module within the software. Place the reaction tubes into the correct positions in your q16 as defined by the software, this may include positioning of empty tubes to ensure that the q16 lid is balanced. The run can then be started.

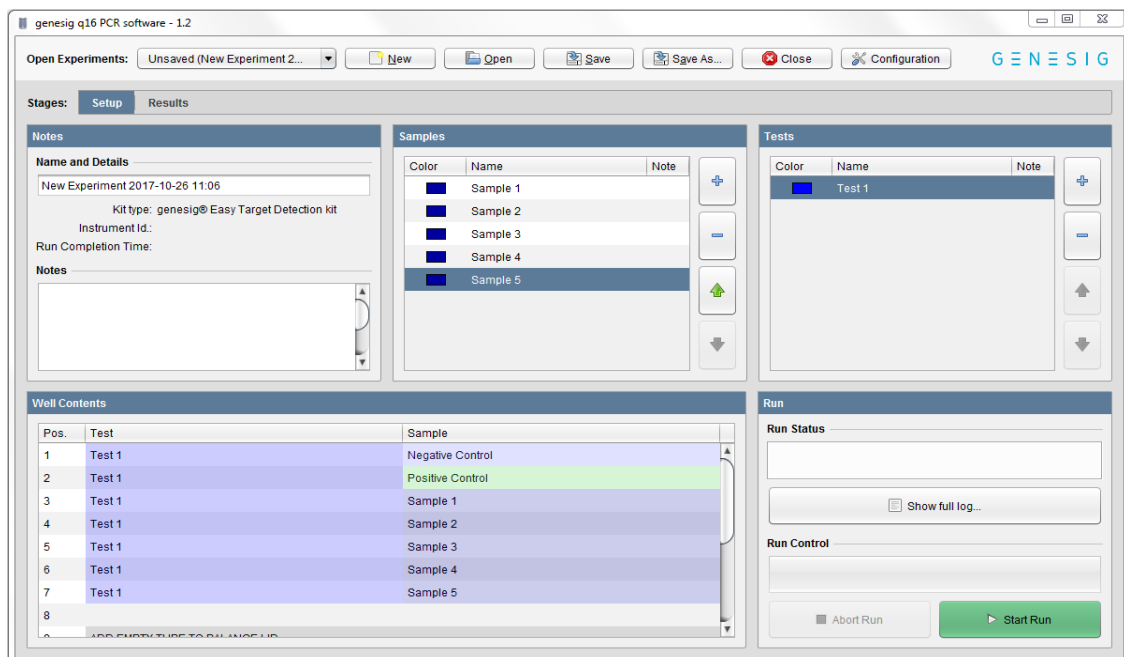

Top tip

- Before loading tubes into the q16, check for bubbles! Flick the bottom of the tubes to remove any bubbles that may have formed during the test setup.
- Apply centrifugal force with a sharp wrist action to ensure all solution is at the bottom of the reaction tube.

What do my results mean?

Analysis of your data is carried out automatically by the genesig® q16 software. The following information is designed to help you fully understand a result or to troubleshoot:

“Positive”

Explanation

Your sample has produced a positive result. Your target of interest is present and you can use the reported quantity.

“Negative”

Explanation

Your sample has produced a negative result. The target is not present in your sample.

“Test contaminated”

Explanation

The Negative Control should be completely free of any DNA/RNA. If you see this error message it means that at some point during the setup, the Negative Control has been contaminated with DNA/RNA and has given a positive signal. This contamination has invalidated the test. The Positive Control and your test samples are both possible sources of contaminating DNA/RNA. The genesig® q16 reaction tubes from previous runs will also contain very high amounts of DNA so it is important that these are carefully disposed of after the run is completed and NEVER OPENED. It may be the case that your kits have become contaminated which will lead to the same problem occurring repeatedly.

Solutions

1. Clean your working area using a commercial DNA remover solution to ensure the area is DNA free at the start of your run and re-run the test.
2. If the problem persists, then the kit has become contaminated, and it will have to be discarded and replaced with a new kit. When you open the new kit, run a simple test to show that changing the kit has solved the problem. Prepare a test which includes only the Positive Control, the Negative Control and one ‘mock sample’. For the ‘mock sample’ add internal control template instead of any sample RNA. The result for the Negative Control and the mock sample should be negative indicating that contamination is no longer present.

Preventive action

An ideal lab set-up has a 'Clean area' where the test reagents are prepared and a 'sample area' where DNA/RNA samples and the Positive Control template are handled. The best workflow involves setting up all the test components (excluding the positive control template) in the clean area and then moving the tests to the sample area for sample and Positive Control addition. If this method is followed, then the kit components are always kept away from possible sources of contamination. For extra security the Negative Control can be completely prepared and sealed in the clean area. All work areas should be decontaminated regularly with DNA remover.

“Sample preparation failed”

Explanation

The test has failed because the quality of the sample was not high enough. The Internal Extraction Control component identifies whether the sample has been prepared correctly and is of suitable quality. This error message means that this quality control test has failed, and the sample quality is not high enough for analysis.

Solutions

1. Check the sample preparation protocol for any user errors then repeat.
2. Poor quality samples can result from overloading the sample preparation protocol with too much starting material. Try reducing the amount of starting material then repeat.
3. Failing to add the Internal Extraction Control RNA to your sample during the sample preparation protocol can also lead to a reported result of “sample preparation failed”. Ensure that this step has not been overlooked or forgotten. If your samples are derived from an archive store or from a process separate from your genesig® Easy extraction kit; you must add 5 µl of Internal Extraction Control RNA template into each 0.5 ml of your sample to make it suitable for use on the q16.

“Positive result, poor quality sample”

Explanation

The test is positive so if you are only interested in obtaining a 'present or absent' answer for your sample then your result is reliable. However, the test contains an Internal Extraction Control component that identifies if the sample is of high quality. This quality control test has failed and the sample is therefore not of high enough quality to accurately calculate the exact copy number of RNA present. If you require quantitative information for your sample then proceed with the solutions mentioned above under “Sample preparation failed”.

“Test failed”

Explanation

The test has failed because the Positive Control has not worked. The Positive Control is present to show that all aspects of the test are working correctly together. When this control test fails, the test as a whole is invalidated. This finding indicates that a problem has occurred in the reaction set-up part of the experiment and has nothing to do with sample preparation.

Solutions

1. Check the entire workflow and test set-up to look for any user errors, then repeat the test e.g., have the right colour pipettes and solutions been used with the correct tubes?
2. Ensure the positive and negative controls are inserted into the correct wells of your q16.
3. A component of the test may have ‘gone off’ due to handling errors, incorrect storage or exceeding the shelf life. When you open a new kit, run a simple test to show that changing the kit has solved the problem. Prepare a test which includes only the Positive Control, the Negative Control and one ‘mock sample’. For the ‘mock sample’ add internal control template instead of any sample RNA. If the Positive Control works, the mock sample will now be called as a negative result.

“Test failed and is contaminated”

Explanation

The Positive Control is indicating test failure, and the Negative Control is indicating test contamination. Please read the “Test Failed” and “Test contamination” sections of this technical support handbook for further explanation.

Solution

For appropriate solutions, read both the “Test failed” and “Test contaminated” sections of this handbook.

Kit storage and stability

This kit is stable for shipping at ambient temperature but should be stored at -20°C upon arrival. Once the lyophilised components have been resuspended, they should not be exposed to temperatures above -20°C for longer than 30 minutes at a time and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances. Primer Design Ltd does not recommend using the kit after the expiry date stated on the pack.

Suitable sample material

This kit can be used with all types of samples from various origins. Please ensure that the extracted nucleic acid samples are suitable in terms of purity, concentration, and RNA integrity.

Dynamic range of test

Under optimal PCR conditions the kit can achieve priming efficiencies between 90-110% and detect less than 100 copies of target template.

Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the USA or the appropriate regulatory authorities in the country of use. During the warranty period Primer Design Ltd genesig® detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations the right to claim under guarantee is expired.

Trademarks

genesig® is a registered trademark of Primer Design Ltd.

oasig® is a trademark of Primer Design Ltd.