

Primerdesign™ Ltd

R00982

# Influenza A virus subtype H5N8 (avian influenza)

**Kit version: 2**

Haemagglutinin gene  
Neuraminidase gene

genesig® Easy Kit

for use on the genesig® q16


50 reactions

GENESIG

Kits by Primerdesign



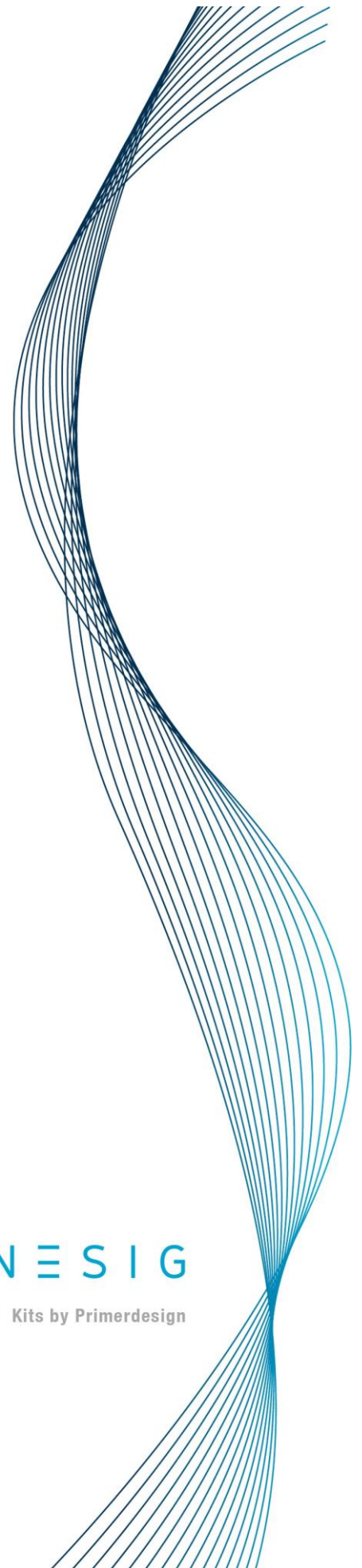

**Specificity of primers and probes last reviewed on:**  
16<sup>th</sup> December 2020

For general laboratory and research use only







# genesig® Easy: at a glance guide





## For each RNA test

Component	Volume	Lab-in-a-box pipette	
H5(H5N8) or N8(H5N8) primer/probe mix	5 µl		
Your RNA sample	5 µl		
oasig™ OneStep Master Mix	10 µl		

## For each positive control

Component	Volume	Lab-in-a-box pipette	
H5(H5N8) or N8(H5N8) primer/probe mix	5 µl		
<u>Positive control template</u>	5 µl		
oasig™ OneStep Master Mix	10 µl		

## For each negative control

Component	Volume	Lab-in-a-box pipette	
H5(H5N8) or N8(H5N8) primer/probe mix	5 µl		
<u>Water</u>	5 µl		
oasig™ OneStep Master Mix	10 µl		

# Kit Contents



- H5(H5N8)-specific primer/probe mix (**BROWN lid**)
  - N8(H5N8)-specific primer/probe mix (**BROWN lid**)
- Once resuspended, the kits should remain at -20°C until ready to use.



- 2x Lyophilised oasig™ OneStep Master Mix




- 2x oasig™ resuspension buffer (**BLUE lid**)



- H5(H5N8) positive control template (**RED lid**)
- N8(H5N8) positive control template (**RED lid**)
- Internal extraction control RNA template (**BLUE lid**)



- RNase/DNase-free water (**WHITE lid**)

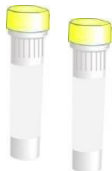

- 2x Template preparation buffer (**YELLOW lid**)



- 2x 54 genesig®q16 reaction tubes

## Reagents and equipment to be supplied by the user

### genesig® q16 instrument

### genesig® Easy Extraction Kit

This kit is designed to work well with all processes that yield high-quality RNA and DNA, but the genesig® Easy extraction method is recommended for ease of use.


### genesig® Lab-In-A-Box

The genesig® Lab-In-A-Box contains all of the pipettes, tips and racks that you will need to use a genesig® Easy kit. Alternatively, if you already have these components and equipment, these can be used instead.

# Step-by-step guide


## 1. Resuspend the test components

Complete this step for each primer/probe mix.

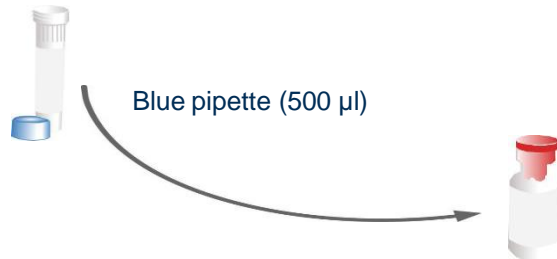

Use the blue pipette to transfer 500 µl\* of the oasisig™ OneStep Master Mix resuspension buffer into the tube of lyophilised oasisig™ OneStep Master Mix and mix well by gently swirling.

\*Transferring 525 µl of the oasisig™ resuspension buffer to your oasisig™ OneStep Master Mix (instead of the 500 µl recommended above) will enable you to take full advantage of the 50 reactions by accounting for volume losses during pipetting. In order to do so with the genesig® Easy fixed volume pipettes use 1x blue, 2x red and 1x grey pipettes to make the total volume. Please be assured that this will not adversely affect the efficiency of the test.

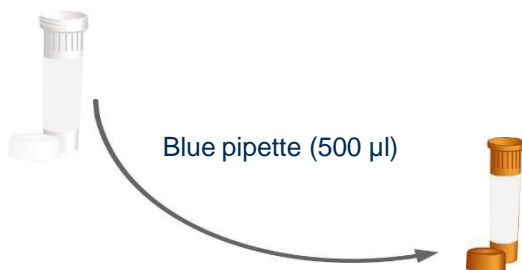


Then use the blue pipette to transfer 500 µl of water into one of the brown tubes labelled H5 or N8 primers/probe. Cap and shake tube to mix. A thorough shake is essential to ensure that all components are resuspended. **Failure to mix well can produce poor kit performance.**

These components are now ready to use. Store them in the freezer from hereon.

### Top tip

- Ensure that the primer/probe mixes are mixed thoroughly before each use by shaking.
- Once resuspended do not expose genesig® Easy kit to temperatures above -20°C for longer than 30 minutes at a time.

## 2. Internal extraction control

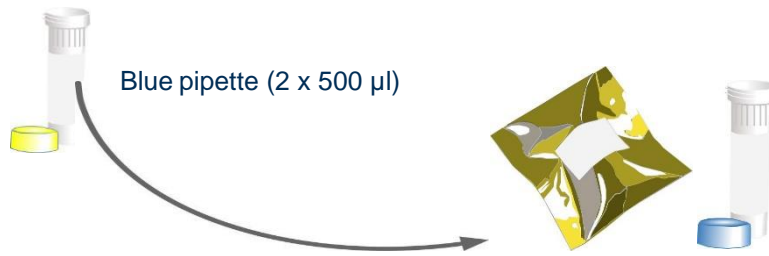

Use the blue pipette to transfer 1000  $\mu\text{l}$  (2 x 500  $\mu\text{l}$ ) of template preparation buffer into the Internal Extraction Control RNA tube. Cap and shake tube to mix.

Your kit contains Internal Extraction Control RNA. This is added to your biological sample at the beginning of the RNA extraction process. It is extracted along with the RNA from your target of interest. The q16 will detect the presence of this Internal Extraction Control RNA at the same time as your target. This is the ideal way to show that your RNA extraction process has been successful.


### **If you are using an alternative extraction kit:**

Use the red pipette to transfer 10  $\mu\text{l}$  of Internal Extraction Control RNA to your sample **after** the lysis buffer has been added, then follow the rest of the extraction protocol.

### **If using samples that have already been extracted:**

Use the grey pipette to transfer 5  $\mu\text{l}$  of Internal Extraction Control RNA to your extracted sample.

## 3. Add primer/probe mix to all reaction tubes

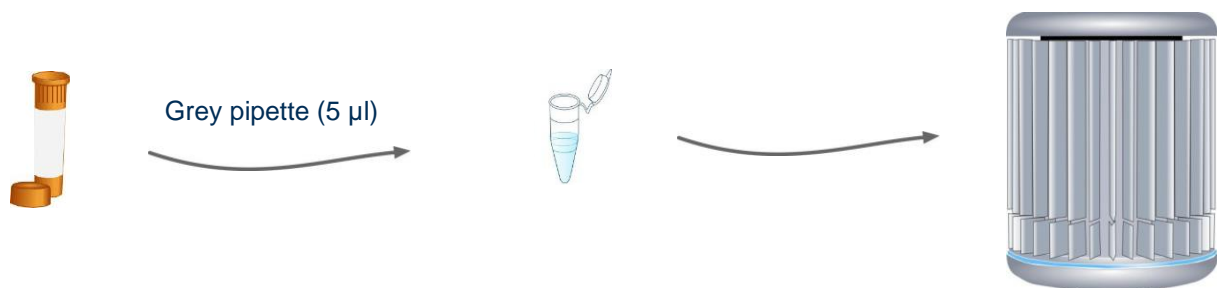


For every reaction run, use the grey pipette to add 5  $\mu\text{l}$  of your H5/N8 primers/probe mix to every tube.

### Top tip

- Always pipette the primer/probe mix directly into the bottom of the tube.
- You can label the tube lids to aid your reaction set-up but avoid labelling tube sides.

## 4. Add Master Mix to all reaction tubes

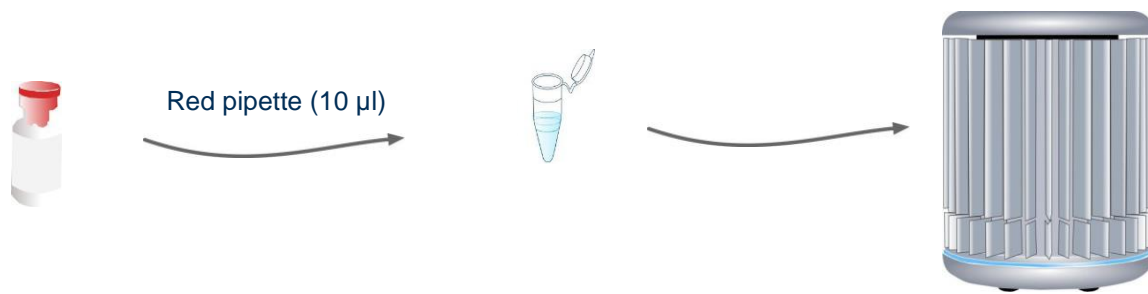


For every reaction run, use the red pipette to add 10 µl of the oasisig™ OneStep Master Mix to the tubes containing primer/probe mix.

Move swiftly to begin your q16 run, as any delay after the oasisig™ OneStep Master Mix has been added can affect the sensitivity of your test.

### Top tip

- Always add the oasisig™ OneStep Master Mix to the side of the tube to reduce the introduction of bubbles.

## 5. Negative control

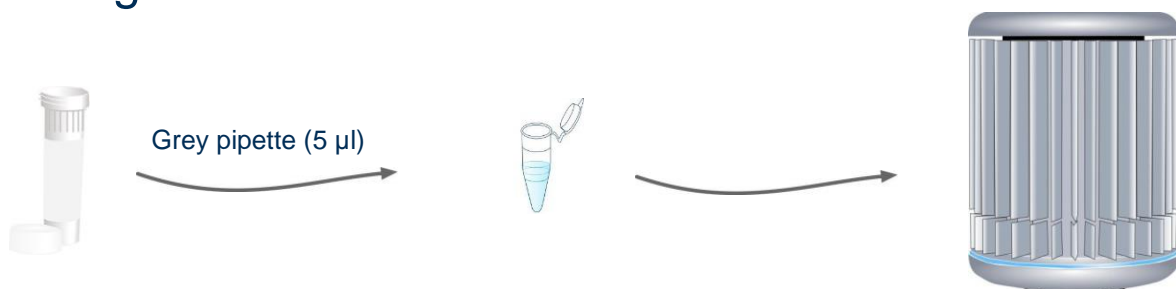

For each test, you will require a negative control. Instead of RNA, water is used. This sample should typically prove negative, thus proving that all of your positive samples are positive.


To create a negative control reaction, simply use the grey pipette to add 5 µl of the water to the required reaction tubes. Close these tubes after adding the water.

Because some genesig® kit targets are common in the environment, you may occasionally see a “late” signal in the negative control. The q16 software will take this into account accordingly.

### Top tip

- Always add the water to the side of the tube to reduce the introduction of bubbles.

## 6. Set up a test

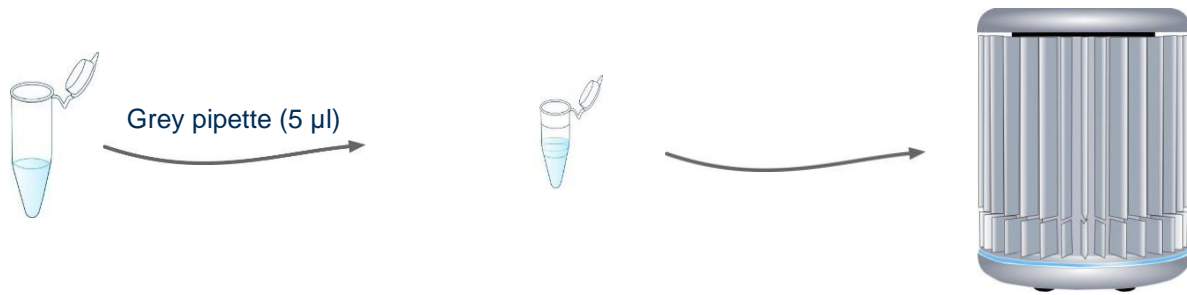


For each sample, you wish to analyse, use the grey pipette to add 5 µl of your RNA sample to the required reaction tubes. Close these tubes after adding the sample. Always change pipette tips between samples.

### Top tip

- Always add the RNA sample to the side of the tube to reduce the introduction of bubbles.

## 7. Positive control

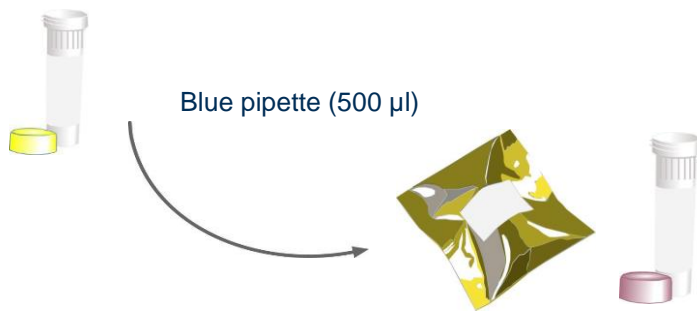


Use the blue pipette to transfer 500 µl of template preparation buffer into the positive control template tube. Cap and shake tube to mix.

Each time you run a test, you will require a positive control. It serves two purposes:

1. It will always test positive to show that everything is working as it should be.
2. The q16 software knows how much RNA is present in the positive control. So, it can automatically compare your sample of interest with the positive control to calculate the amount of target RNA in your sample.

To create a positive control reaction, simply use 5µl of the positive control instead of your RNA sample.

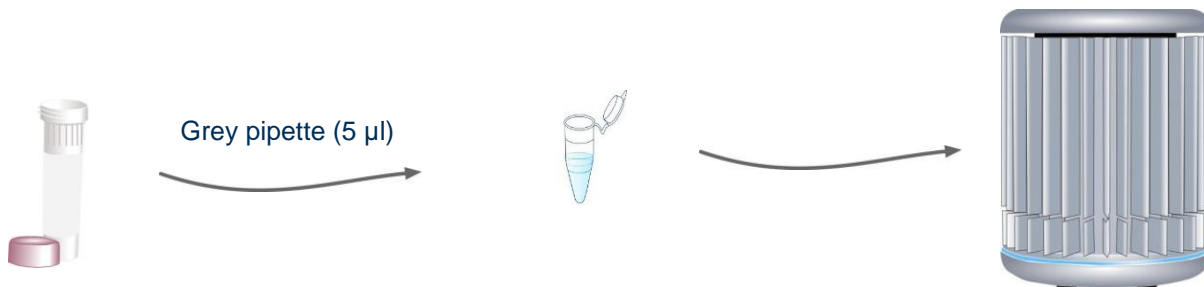


Take great care when setting up your positive control. The positive control template has the potential to give you a false positive signal in your other samples. Set up the positive controls last after all other sample tubes are closed. Always change pipette tips between samples. You may even choose to set up positive controls in a separate room.

### Top tip

- Always add the positive control template to the side of the tube to reduce the introduction of bubbles.

## 8. Running the test

Place the tubes into the correct positions in your q16 as defined by the software; this may include positioning empty tubes to ensure that the q16 lid is balanced. The run can then be started.

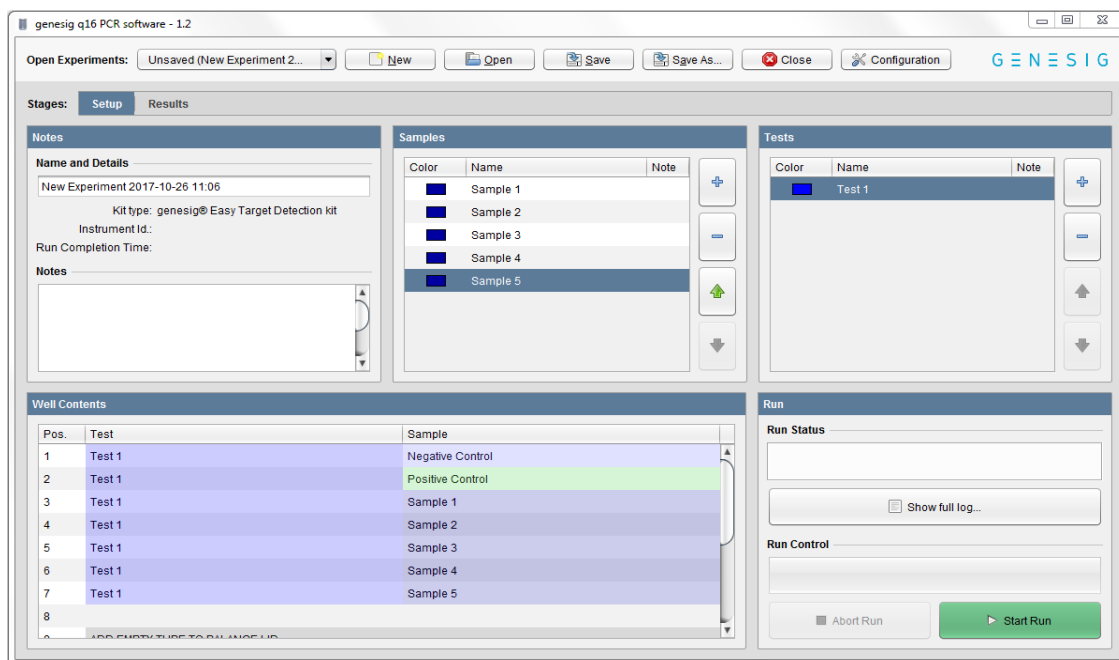

### Top tip

- Before loading tubes into the q16, check for bubbles! Flick the bottom of the tubes to remove any bubbles that may have formed during the test set-up.
- Apply centrifugal force with a sharp wrist action to ensure all solution is at the bottom of the reaction tube.
- When repeating a test, you can use a previous file as a template by clicking 'open' then selecting File name > Files of Type > Experiment file as template.



# What do my results mean?

Analysis of your data is carried out automatically by the q16. The following information is designed to help you fully understand a result or to troubleshoot:

## “Positive”

### Explanation

Your sample has produced a positive result. Your target of interest is present, and you can use the reported quantity.

## “Negative”

### Explanation

Your sample has produced a negative result. The target is not present in your sample.

## “Test contaminated”

### Explanation

The Negative Control should be completely free of any DNA/RNA. If you see this error message, it means that at some point during the set-up, the Negative Control has been contaminated with DNA/RNA and has given a positive signal. This contamination has invalidated the test. The Positive Control and your test samples are both possible sources of contaminating DNA/RNA. The genesig® q16 reaction tubes from previous runs will also contain very high amounts of DNA, so it is important that these are carefully disposed of after the run is completed and NEVER OPENED. It may be the case that your kits have become contaminated, which will lead to the same problem occurring repeatedly.

### Solutions

1. Clean your working area using a commercial DNA remover solution to ensure the area is DNA-free at the start of your run and re-run the test.
2. If the problem persists, then the kit has become contaminated, and it will have to be discarded and replaced with a new kit. When you open the new kit, run a simple test to show that changing the kit has solved the problem. Prepare a test that includes only the Positive Control, the Negative Control and one ‘mock sample’. For the ‘mock sample’ add water instead of any sample RNA. The result for the Negative Control and the mock sample should be negative, indicating that contamination is no longer present.

### Preventive action

An ideal lab set-up has a ‘Clean area’ where the test reagents are prepared and a ‘sample area’ where DNA/RNA samples and the Positive Control template are handled. The best workflow involves setting up all the test components (excluding the positive control template) in the clean area and then moving the tests to the sample area for the sample and Positive Control addition.

## “Sample preparation failed”

### Explanation

The test has failed because the quality of the sample was not high enough. The Internal Extraction Control component identifies whether the sample has been prepared correctly and is of suitable quality. This error message means that this quality control test has failed, and the sample quality is not high enough for analysis.

### Solutions

1. Check the sample preparation protocol for any user errors, then repeat.
2. Poor quality samples can result from overloading the sample preparation protocol with too much starting material. Try reducing the amount of starting material then repeat.
3. Failing to add the Internal Extraction Control RNA to your sample during the sample preparation protocol can also lead to a reported result of “sample preparation failed”. Ensure that this step has not been overlooked or forgotten. If your samples are derived from an archive store or from a process separate from your genesig® Easy extraction kit; you must add 5µl of Internal Extraction Control RNA into each 0.5 ml of your sample to make it suitable for use on the q16.

## “Positive result, poor quality sample”

### Explanation

The test is positive, so if you are only interested in obtaining a ‘present or absent’ answer for your sample, then your result is reliable. However, the test contains an Internal Extraction Control component that identifies if the sample is of high quality. This quality control test has failed, therefore the sample is not of high enough quality to accurately calculate the exact copy number of RNA present. If you require quantitative information for your sample, then proceed with the solutions below.

## “Test failed”

### Explanation

The test has failed because the Positive Control has not worked. The Positive Control is present to show that all aspects of the test are working correctly together. When this control test fails, the test as a whole is invalidated. This finding indicates that a problem has occurred in the reaction set-up part of the experiment and has nothing to do with sample preparation.

### Solutions

1. Check the entire workflow and test set-up to look for any user errors, then repeat the test e.g., have the right colour pipettes and solutions been used with the correct tubes?
2. Ensure the positive and negative controls are inserted into the correct wells of your q16.
3. A component of the test may have ‘gone off’ due to handling errors, incorrect storage or exceeding the shelf life. When you open a new kit, run a simple test to show that changing the kit has solved the problem. Prepare a test that includes only the Positive Control, the Negative Control and one ‘mock sample’. For the ‘mock sample’ add internal control template instead of any sample RNA. If the Positive Control works, the mock sample will now be called as a negative result.

## “Test failed and is contaminated”

### Explanation

The Positive Control is indicating test failure, and the Negative Control is indicating test contamination. Please read the “Test Failed” and “Test contamination” sections of this technical support handbook for a further explanation.

### Solution

1. For appropriate solutions, read both the “Test failed” and “Test contaminated” sections of this handbook.

## Disposal Considerations

Dispose of tested samples according to any local, national, or regional regulations.

Product: Offer surplus and non-recyclable solutions to a licensed disposal company.

Contaminated packaging: Dispose of as unused product.

# Introduction to Influenza A virus subtype H5N8

Influenza type A viruses are 80–120 nanometers in diameter and usually roughly spherical, made up of a viral envelope containing two main types of proteins, wrapped around a central core. The two large proteins found on the outside of viral particles are haemagglutinin (HA) and neuraminidase (NA). HA is a protein that mediates binding of the virion to target cells and entry of the viral genome into the target cell, while NA is involved in the release of progeny virions from infected cells. Influenza type A viruses are categorized into subtypes based on the type of these two proteins on the surface of the viral envelope. The central core of a virion contains the viral genome and other viral proteins that package and protect the genetic material. The entire Influenza A virus genome is ~13,588 bases long and is contained on 8 RNA segments that code for 11 proteins.

Avian Influenza ('avian flu' or 'bird flu') strains all belong to the Influenza A virus which are generally adapted to birds but not exclusive to them. Of particular concern to humans are ones which are infectious to both humans and birds. Specific strains such as H1N1 have been the subject of much media concern and speculation over pandemics and its widespread transmission globally. Only two influenza A virus subtypes A(H1N1)pdm09, and A(H3N2), are currently circulating among people..

Since 2016, H5N8 has been spreading rapidly through migration of wild birds in Europe and Asia, causing deaths in the wild bird population and domestic poultry. An outbreak of H5N8 in Russia was confirmed by WHO in 2020, with the virus also detected in poultry or wild birds in Bulgaria, the Czech Republic, Egypt, Germany, Hungary, Iraq, Japan, Kazakhstan, the Netherlands, Poland, Romania, and the United Kingdom. There is a low possibility of human infection with H5N8, however, it should not be excluded as infections in seven humans were reported by Russia in 2021, the first confirmed cases of H5N8 in humans as a result of the 2020 outbreak.

# Specificity

The Primerdesign™ genesig® Kit for Influenza A virus subtype H5N8 (avian influenza) is designed for the in vitro quantification of H5N8 genomes. The kit is designed to have a broad detection profile. Specifically, the primers will detect over 95% of the sequences on the GISAID EpiFlu database, identified as H5N8 from avian origin and collected within the three years to the time of the most recent review. Due to the sequence evolution observed in the genomes of RNA viruses, a period of three years is considered an adequate time period for analysis.

The H5 primer and probe set is designed to detect all relevant H5 sequences from the H5N8 lineage but may cross react with H5 sequences from other subtypes.

The N8 primer and probe set is designed to detect all relevant N8 sequences from the H5N8 lineage but may cross react with N8 sequences from other subtypes.

The dynamics of genetic variation mean that new sequence information may become available after the most recent review. Primerdesign periodically reviews the detection profiles of our kits and when required releases new versions.

If you require further information or have a specific question about the detection profile of this kit then please send an e-mail to [techsupport@primerdesign.co.uk](mailto:techsupport@primerdesign.co.uk) and our team will answer your question.

## Kit storage and stability

This lyophilised kit is stable at room temperature but should be stored at -20°C on arrival. Once the lyophilised components have been resuspended, they should not be exposed to temperatures above -20°C for longer than 30 minutes at a time and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances.

Primerdesign does not recommend using the kit after the expiry date stated on the pack.

## Suitable sample material

All kinds of sample material suited for PCR amplification can be used. Please ensure the samples are suitable in terms of purity, concentration, and RNA/DNA integrity.

## Dynamic range of test

Under optimal PCR conditions the kit can achieve priming efficiencies between 90-110% and detect less than 100 copies of target template.

## Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the USA or the appropriate regulatory authorities in the country of use. During the warranty period, Primerdesign™ genesig® detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations, the right to claim under guarantee is expired. PCR is a proprietary technology covered by several US and foreign patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation in certain fields. Depending on your specific application, you may need a license from Roche or PE to practice PCR. Additional information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Roche Molecular Systems, 1145 Atlantic Avenue, Alameda, CA 94501 or Applied Biosystems business group of the Applied Biosystems Corporation, 850 Lincoln Centre Drive, Foster City, CA 94404. In addition, the 5' nuclease assay and other homogeneous amplification methods used in connection with the PCR process may be covered by US Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc, and by US Patent 5,538,848, owned by The Perkin-Elmer Corporation.

## Trademarks

Primerdesign™ is a trademark of Primerdesign Ltd.  
genesig® is a registered trademark of Primerdesign Ltd.  
oasig™ is a trademark of Primerdesign Ltd.