

Primerdesign™ Ltd

Bovine adenovirus 5/6/8

Hexon gene

genesig® Standard Kit

150 tests

G E N E S I G

Kits by Primerdesign

For general laboratory and research use only

Introduction to Bovine adenovirus 5/6/8

Bovine Adenovirus (BAdV) belongs to the Adenoviridae family and is found worldwide. Currently there are 10 recognised serotypes of adenovirus in cattle, as an adenovirus BAdV is approx. 34 kbp in length with a G+C content of ~54%.

Bovine adenoviruses (BAdVs) belong to the Mastadenovirus genus (which includes adenoviruses of mammals) of the family Adenoviridae and are involved in respiratory and enteric infections of calves. Bovine Adenovirus can infect cattle and zebu of any breed, sex or age, however, younger animals are more at risk as their maternal antibody begins to wane, around the age of 2 weeks to 4 months old. Experimentally, the virus has also been shown to infect goats and sheep.

Infection usually results in disease of the gastrointestinal or respiratory tract. However, infection may not always result in disease as the virus can be isolated in healthy cattle.

Once infected, the cattle shed the virus for approximately 10 days in the respiratory secretions or faeces -- some cattle may become persistently infected, resulting in excretion of the virus for much longer.

Symptoms of the gastrointestinal tract include diarrhoea with mucus, melena or haematochezia, anorexia or reduced appetite, abdominal distension and possible dysphagia. Respiratory signs include coughing (sometimes with blood), serous nasal discharge, dyspnoea and tachypnoea. This may progress to bronchopneumonia if secondary bacterial infection is present.

Generalised clinical symptoms include pyrexia, weight loss, sudden death, lymphadenopathy and generalised weakness. The animal will also be dull, depressed and lethargic.

Clinical signs and physical examination may add Bovine Adenovirus to the differential list, but definitive diagnosis cannot be achieved without clinical testing. Because adenoviruses can be isolated from apparently healthy cattle, isolation of adenoviruses from clinically sick calves does not necessarily mean that the isolated adenovirus type plays an aetiological role in the clinical disease under consideration.

Because of the number of types of adenoviruses infecting cattle, virus isolation is necessary to definitively identify the virus. Virus can be isolated from nasal secretions, tracheal fluids, intestinal contents and tissue homogenates and used in nucleic acid testing via Real-Time PCR.

When economically feasible, cattle can be treated to provide relief from clinical signs associated with adenoviral infection.

Because secondary bacterial infections are common in cattle, antibiotics are often used as part of the treatment.

The main control measure is to ensure adequate colostrum at birth as passive transfer provides immunity to calves.

Both modified live and inactivated adenovirus vaccines have been developed and evaluated for use in cattle and they should be administered when maternal antibodies have waned, but 2 to 3 weeks before calves from different places are assembled under stressful conditions. Such vaccines are available in Europe, though not in the UK and Japan.

Specificity

The Primerdesign genesig Kit for Bovine adenovirus 5/6/8 (BAdV-5/6/8) genomes is designed for the in vitro quantification of BAdV-5/6/8 genomes. The kit is designed to have a broad detection profile. Specifically, the primers represent 100% homology with over 95% of the NCBI database reference sequences available at the time of design.

The dynamics of genetic variation means that new sequence information may become available after the initial design. Primerdesign periodically reviews the detection profiles of our kits and when required releases new versions.

If you require further information, or have a specific question about the detection profile of this kit then please send an e.mail to enquiry@primerdesign.co.uk and our bioinformatics team will answer your question.

Kit contents

- **BAdV-5/6/8 specific primer/probe mix (150 reactions BROWN)**
FAM labelled
- **BAdV-5/6/8 positive control template (for Standard curve RED)**
- **RNase/DNase free water (WHITE)**
for resuspension of primer/probe mixes
- **Template preparation buffer (YELLOW)**
for resuspension of positive control template and standard curve preparation

Reagents and equipment to be supplied by the user

Real-time PCR Instrument

Extraction kit

This kit is recommended for use with genesig Easy DNA/RNA extraction kit. However, it is designed to work well with all processes that yield high quality RNA and DNA with minimal PCR inhibitors.

oasig™ lyophilised or Precision® PLUS 2X qPCR Master Mix

This kit is intended for use with oasig or PrecisionPLUS 2X qPCR Master Mix.

Pipettors and Tips

Vortex and centrifuge

Thin walled 1.5 ml PCR reaction tubes

Kit storage and stability

This kit is stable at room temperature but should be stored at -20°C on arrival. Once the lyophilised components have been resuspended they should not be exposed to temperatures above -20°C for longer than 30 minutes at a time and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances.

If a standard curve dilution series is prepared this can be stored frozen for an extended period. If you see any degradation in this serial dilution a fresh standard curve can be prepared from the positive control.

Primerdesign does not recommend using the kit after the expiry date stated on the pack.

Suitable sample material

All kinds of sample material suited for PCR amplification can be used. Please ensure the samples are suitable in terms of purity, concentration, and DNA integrity. Always run at least one negative control with the samples. To prepare a negative-control, replace the template DNA sample with RNase/DNase free water.

Dynamic range of test

Under optimal PCR conditions genesig BAdV-5/6/8 detection kits have very high priming efficiencies of >95% and can detect less than 100 copies of target template.

Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the USA or the appropriate regulatory authorities in the country of use. During the warranty period Primerdesign genesig detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations the right to claim under guarantee is expired. PCR is a proprietary technology covered by several US and foreign patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation in certain fields. Depending on your specific application you may need a license from Roche or PE to practice PCR. Additional information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Roche Molecular Systems, 1145 Atlantic Avenue, Alameda, CA 94501 or Applied Biosystems business group of the Applied Biosystems Corporation, 850 Lincoln Centre Drive, Foster City, CA 94404. In addition, the 5' nuclease assay and other homogeneous amplification methods used in connection with the PCR process may be covered by U.S. Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc., and by U.S. Patent 5,538,848, owned by The Perkin-Elmer Corporation.

Trademarks

Primerdesign™ is a trademark of Primerdesign Ltd.

genesig® is a registered trademark of Primerdesign Ltd.

The PCR process is covered by US Patents 4,683,195, and 4,683,202 and foreign equivalents owned by Hoffmann-La Roche AG. BI, ABI PRISM®, GeneAmp® and MicroAmp® are registered trademarks of the Applied Biosystems (Applied Biosystems Corporation). BIOMEK® is a registered trademark of Beckman Instruments, Inc.; iCycler™ is a registered trademark of Bio-Rad Laboratories, Rotor-Gene is a trademark of Corbett Research. LightCycler™ is a registered trademark of the Idaho Technology Inc. GeneAmp®, TaqMan® and AmpliTaqGold® are registered trademarks of Roche Molecular Systems, Inc., The purchase of the Primerdesign reagents cannot be construed as an authorization or implicit license to practice PCR under any patents held by Hoffmann-La Roche Inc.

Principles of the test

Real-time PCR

A BAdV-5/6/8 specific primer and probe mix is provided and this can be detected through the FAM channel.

The primer and probe mix provided exploits the so-called TaqMan® principle. During PCR amplification, forward and reverse primers hybridize to the BAdV-5/6/8 DNA. A fluorogenic probe is included in the same reaction mixture which consists of a DNA probe labeled with a 5`-dye and a 3`-quencher. During PCR amplification, the probe is cleaved and the reporter dye and quencher are separated. The resulting increase in fluorescence can be detected on a range of qPCR platforms.

Positive control

For copy number determination and as a positive control for the PCR set up, the kit contains a positive control template. This can be used to generate a standard curve of BAdV-5/6/8 copy number / Cq value. Alternatively the positive control can be used at a single dilution where full quantitative analysis of the samples is not required. Each time the kit is used, at least one positive control reaction must be included in the run. A positive result indicates that the primers and probes for detecting the target BAdV-5/6/8 gene worked properly in that particular experimental scenario. If a negative result is obtained the test results are invalid and must be repeated. Care should be taken to ensure that the positive control does not contaminate any other kit component which would lead to false-positive results. This can be achieved by handling this component in a Post PCR environment. Care should also be taken to avoid cross-contamination of other samples when adding the positive control to the run. This can be avoided by sealing all other samples and negative controls before pipetting the positive control into the positive control well.

Negative control

To validate any positive findings a negative control reaction should be included every time the kit is used. For this reaction the RNase/DNase free water should be used instead of template. A negative result indicates that the reagents have not become contaminated while setting up the run.

Resuspension protocol

To minimize the risk of contamination with foreign DNA, we recommend that all pipetting be performed in a PCR clean environment. Ideally this would be a designated PCR lab or PCR cabinet. Filter tips are recommended for all pipetting steps.

1. Pulse-spin each tube in a centrifuge before opening.

This will ensure lyophilised primer and probe mix is in the base of the tube and is not spilt upon opening the tube.

2. Resuspend the kit components in the RNase/DNase free water supplied, according to the table below.

To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in water	Volume
Pre-PCR pack	
BAdV-5/6/8 primer/probe mix (BROWN)	165 µl

3. Resuspend the positive control template in the template preparation buffer supplied, according to the table below:

To ensure complete resuspension, vortex the tube thoroughly.

Component - resuspend in template preparation buffer	Volume
Post-PCR heat-sealed foil	
BAdV-5/6/8 Positive Control Template (RED) *	500 µl

* This component contains high copy number template and is a VERY significant contamination risk. It must be opened and handled in a separate laboratory environment, away from the other components.

qPCR detection protocol

1. **For each DNA sample prepare a reaction mix according to the table below:**
Include sufficient reactions for positive and negative controls.

Component	Volume
oasig or PrecisionPLUS 2X qPCR Master Mix	10 μ l
BAdV-5/6/8 primer/probe mix (BROWN)	1 μ l
RNase/DNase free water (WHITE)	4 μ l
Final Volume	15 μl

2. **Pipette 15 μ l of this mix into each well according to your qPCR experimental plate set up.**
3. **Prepare DNA templates for each of your samples.**
4. **Pipette 5 μ l of DNA template into each well, according to your experimental plate set up.**
For negative control wells use 5 μ l of RNase/DNase free water. The final volume in each well is 20 μ l.
5. **If a standard curve is included for quantitative analysis, prepare a reaction mix according to the table below:**

Component	Volume
oasig or PrecisionPLUS 2X qPCR Master Mix	10 μ l
BAdV-5/6/8 primer/probe mix (BROWN)	1 μ l
RNase/DNase free water (WHITE)	4 μ l
Final Volume	15 μl

6. Preparation of a standard curve dilution series.

- 1) Pipette 90µl of template preparation buffer into 5 tubes and label 2-6
- 2) Pipette 10µl of Positive Control Template (RED) into tube 2
- 3) Vortex thoroughly
- 4) Change pipette tip and pipette 10µl from tube 2 into tube 3
- 5) Vortex thoroughly

Repeat steps 4 and 5 to complete the dilution series

Standard Curve	Copy Number
Tube 1 Positive control (RED)	2×10^5 per µl
Tube 2	2×10^4 per µl
Tube 3	2×10^3 per µl
Tube 4	2×10^2 per µl
Tube 5	20 per µl
Tube 6	2 per µl

7. Pipette 5µl of standard template into each well for the standard curve according to your experimental plate set up.
The final volume in each well is 20µl.

qPCR amplification protocol

Amplification conditions using oasis or PrecisionPLUS 2X qPCR Master Mix.

	Step	Time	Temp
Cycling x50	Enzyme activation	2 min	95 °C
	Denaturation	10 s	95 °C
	DATA COLLECTION *	60 s	60 °C


* Fluorogenic data should be collected during this step through the FAM channel

Interpretation of results


Target	Positive control	Negative control	Interpretation
+	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
-	+	-	NEGATIVE RESULT
+ / -	+	≤ 35	EXPERIMENT FAILED due to test contamination
+ / -	+	> 35	*
+ / -	-	+ / -	EXPERIMENT FAILED

Positive control template (RED) is expected to amplify between Cq 16 and 23. Failure to satisfy this quality control criterion is a strong indication that the experiment has been compromised

*Where the test sample is positive and the negative control is positive with a Cq > 35 , the sample must be reinterpreted based on the relative signal strength of the two results:

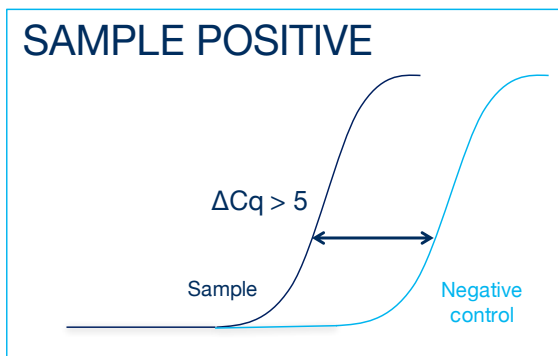

If the sample amplifies > 5 Cq earlier than the negative control then the sample should be reinterpreted (via the table above) with the negative control verified as negative.

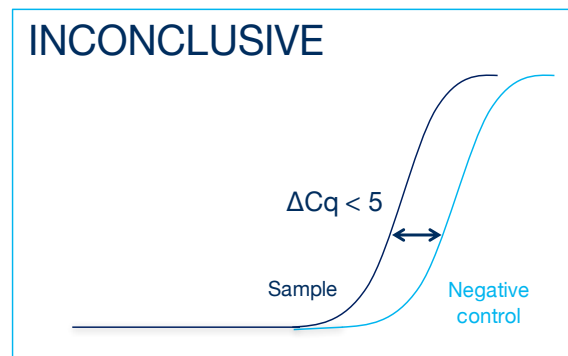

If the sample amplifies < 5 Cq earlier than the negative control then the positive sample result is invalidated and the result should be determined inconclusive due to test contamination. The test for this sample should be repeated.