

Primerdesign™ Ltd

Mumps Virus

small hydrophobic protein (SH) gene

genesig® Advanced Kit

150 tests

GENESIG

Kits by Primerdesign

For general laboratory and research use only

Introduction to Mumps Virus

The mumps virus is the causative agent of the disease mumps. Prior to the introduction of routine vaccination, it was once a common childhood disease worldwide and still is in the third world. Mumps virus belongs to the genus Rubulavirus in the family Paramyxoviridae and is seen to have a roughly spherical, enveloped morphology of about 200 nm in diameter. It contains a non-segmented, linear, single-stranded, negative-sense RNA genome of approximately 15.4 kb that is coated by nucleocapsid proteins. The genome contains open reading frames for a nucleoprotein (N), phosphoprotein (P), matrix protein (M), fusion protein (F), small hydrophobic protein (SH), haemagglutinin–neuraminidase protein (HN) and polymerase protein (L). In addition, the P gene synthesizes the V and I proteins through an RNA-editing mechanism.

Transmission is mainly through either the aerosol route or direct contact but can also be spread by contact with contaminated objects. The virus enters the body through the pharynx or the conjunctiva, the nucleocapsid being released into the host cells through membrane fusion after receptor binding to sialic acid located on the host cells. As such, the virus infects most cell types due to the ubiquitous nature of the receptor on mammalian cells. The virus localises in the throat but disseminates through the bloodstream, allowing it to infect various organs, namely, salivary glands, testes, ovaries, pancreas, and the brain. Viral replication occurs in the ductal cells of the glands, inducing large amounts of complementary RNA in the infected cells which is protected from degradation by nucleoprotein. The virus also protects itself from host defenses through the activity of V protein, which acts on the regulatory RNA helicase MDA5 to circumvent interferon- α and interferon- γ signaling, and a small hydrophobic (SH) protein that interferes with tumor necrosis factor alpha-mediated apoptosis.

The incubation period for mumps virus is usually 12-24 days. Infection can be asymptomatic in 20%-30% of cases. Classical symptoms include fever, headache, face pain, malaise, unilateral or bilateral swelling of the parotid glands and jaw and orchitis in males. The central nervous system is affected in approximately 50% of cases and more serious complications include deafness.

Specificity

The Primerdesign genesig Kit for Mumps Virus (Mumps) genomes is designed for the in vitro quantification of Mumps genomes. The kit is designed to have a broad detection profile. Specifically, the primers represent 100% homology with over 95% of the NCBI database reference sequences available at the time of design.

The dynamics of genetic variation means that new sequence information may become available after the initial design. Primerdesign periodically reviews the detection profiles of our kits and when required releases new versions.

Our kit for Mumps Virus has been designed for the specific and exclusive in vitro quantification of the Mumps Virus genome. The SH target gene has previously been shown to be a good marker for the detection of the Mumps Virus (Boddicker JD et. al., 2007).

The primers and probe sequences in this kit have 100% homology with over 95% of reference sequences in the NCBI database based on a comprehensive bioinformatics analysis.

If you require further information, or have a specific question about the detection profile of this kit then please send an e.mail to enquiry@primerdesign.co.uk and our bioinformatics team will answer your question.

Kit contents

- Mumps specific primer/probe mix (150 reactions **BROWN**)
FAM labelled
- Mumps positive control template (for Standard curve **RED**)
- Internal extraction control primer/probe mix (150 reactions **BROWN**)
VIC labelled as standard
- Internal extraction control RNA (150 reactions **BLUE**)
- Endogenous control primer/probe mix (150 reactions **BROWN**)
FAM labelled
- RNase/DNase free water (**WHITE**)
for resuspension of primer/probe mixes
- Template preparation buffer (**YELLOW**)
for resuspension of internal control template, positive control template and standard curve preparation

Reagents and equipment to be supplied by the user

Real-time PCR Instrument

Extraction kit

This kit is recommended for use with genesig Easy DNA/RNA Extraction kit. However, it is designed to work well with all processes that yield high quality RNA and DNA with minimal PCR inhibitors.

oasig™ lyophilised OneStep or Precision®PLUS OneStep 2X RT-qPCR Master Mix

Contains complete OneStep RT-qPCR master mix

Pipettors and Tips

Vortex and centrifuge

Thin walled 1.5 ml PCR reaction tubes

Kit storage and stability

This kit is stable at room temperature but should be stored at -20°C on arrival. Once the lyophilised components have been resuspended they should not be exposed to temperatures above -20°C for longer than 30 minutes at a time and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances.

If a standard curve dilution series is prepared this can be stored frozen for an extended period. If you see any degradation in this serial dilution a fresh standard curve can be prepared from the positive control.

Primerdesign does not recommend using the kit after the expiry date stated on the pack.

Suitable sample material

All kinds of sample material suited for PCR amplification can be used. Please ensure the samples are suitable in terms of purity, concentration, and RNA/DNA integrity (An internal PCR control is supplied to test for non specific PCR inhibitors). Always run at least one negative control with the samples. To prepare a negative-control, replace the template RNA sample with RNase/DNase free water.

Dynamic range of test

Under optimal PCR conditions genesig Mumps detection kits have very high priming efficiencies of >95% and can detect less than 100 copies of target template.

Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the USA or the appropriate regulatory authorities in the country of use. During the warranty period Primerdesign genesig detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations the right to claim under guarantee is expired. PCR is a proprietary technology covered by several US and foreign patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation in certain fields. Depending on your specific application you may need a license from Roche or PE to practice PCR. Additional information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Roche Molecular Systems, 1145 Atlantic Avenue, Alameda, CA 94501 or Applied Biosystems business group of the Applied Biosystems Corporation, 850 Lincoln Centre Drive, Foster City, CA 94404. In addition, the 5' nuclease assay and other homogeneous amplification methods used in connection with the PCR process may be covered by U.S. Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc. and by U.S. Patent 5,538,848, owned by The Perkin-Elmer Corporation.

Trademarks

Primerdesign™ is a trademark of Primerdesign Ltd.

genesig® is a registered trademark of Primerdesign Ltd.

The PCR process is covered by US Patents 4,683,195, and 4,683,202 and foreign equivalents owned by Hoffmann-La Roche AG. BI, ABI PRISM®, GeneAmp® and MicroAmp® are registered trademarks of the Applied Biosystems (Applied Biosystems Corporation). BIOMEK® is a registered trademark of Beckman Instruments, Inc.; iCycler™ is a registered trademark of Bio-Rad Laboratories, Rotor-Gene is a trademark of Corbett Research. LightCycler™ is a registered trademark of the Idaho Technology Inc. GeneAmp®, TaqMan® and AmpliTaqGold® are registered trademarks of Roche Molecular Systems, Inc., The purchase of the Primerdesign™ reagents cannot be construed as an authorization or implicit license to practice PCR under any patents held by Hoffmann-La Roche Inc.

Principles of the test

Real-time PCR

A Mumps specific primer and probe mix is provided and this can be detected through the FAM channel.

The primer and probe mix provided exploits the so-called TaqMan® principle. During PCR amplification, forward and reverse primers hybridize to the Mumps cDNA. A fluorogenic probe is included in the same reaction mixture which consists of a DNA probe labeled with a 5`-dye and a 3`-quencher. During PCR amplification, the probe is cleaved and the reporter dye and quencher are separated. The resulting increase in fluorescence can be detected on a range of qPCR platforms.

Positive control

For copy number determination and as a positive control for the PCR set up, the kit contains a positive control template.

This can be used to generate a standard curve of Mumps copy number / Cq value. Alternatively the positive control can be used at a single dilution where full quantitative analysis of the samples is not required. Each time the kit is used, at least one positive control reaction must be included in the run. A positive result indicates that the primers and probes for detecting the target Mumps gene worked properly in that particular experimental scenario. If a negative result is obtained the test results are invalid and must be repeated. Care should be taken to ensure that the positive control does not contaminate any other kit component which would lead to false-positive results. This can be achieved by handling this component in a Post PCR environment. Care should also be taken to avoid cross-contamination of other samples when adding the positive control to the run. This can be avoided by sealing all other samples and negative controls before pipetting the positive control into the positive control well.

Negative control

To validate any positive findings a negative control reaction should be included every time the kit is used. For this reaction the RNase/DNase free water should be used instead of template. A negative result indicates that the reagents have not become contaminated while setting up the run.

Internal RNA extraction control

When performing RNA extraction, it is often advantageous to have an exogenous source of RNA template that is spiked into the lysis buffer. This control RNA is then co-purified with the sample RNA and can be detected as a positive control for the extraction process. Successful co-purification and qPCR for the control RNA also indicates that PCR inhibitors are not present at a high concentration.

A separate qPCR primer/probe mix are supplied with this kit to detect the exogenous RNA using qPCR. The PCR primers are present at PCR limiting concentrations which allows multiplexing with the target sequence primers. Amplification of the control cDNA does not interfere with detection of the Mumps target cDNA even when present at low copy number. The Internal control is detected through the VIC channel and gives a Cq value of 28+/-3 depending on the level of sample dilution.

Endogenous control

To confirm extraction of a valid biological template, a primer and probe mix is included to detect an endogenous gene. Detection of the endogenous control is through the FAM channel and it is NOT therefore possible to perform a multiplex with the Mumps primers. A poor endogenous control signal may indicate that the sample did not contain sufficient biological material.

Resuspension protocol

To minimize the risk of contamination with foreign DNA, we recommend that all pipetting be performed in a PCR clean environment. Ideally this would be a designated PCR lab or PCR cabinet. Filter tips are recommended for all pipetting steps.

1. Pulse-spin each tube in a centrifuge before opening.

This will ensure lyophilised primer and probe mix is in the base of the tube and is not spilt upon opening the tube.

2. Resuspend the primer/probe mixes in the RNase/DNase free water supplied, according to the table below:

To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in water	Volume
Pre-PCR pack	
Mumps primer/probe mix (BROWN)	165 µl
Internal extraction control primer/probe mix (BROWN)	165 µl
Endogenous control primer/probe mix (BROWN)	165 µl

3. Resuspend the internal control template and positive control template in the template preparation buffer supplied, according to the table below:

To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in template preparation buffer	Volume
Pre-PCR heat-sealed foil	
Internal extraction control RNA (BLUE)	600 µl
Post-PCR heat-sealed foil	
Mumps Positive Control Template (RED) *	500 µl

* This component contains high copy number template and is a VERY significant contamination risk. It must be opened and handled in a separate laboratory environment, away from the other components.

RNA extraction

The internal extraction control RNA can be added either to the RNA lysis/extraction buffer or to the RNA sample once it has been resuspended in lysis buffer.

DO NOT add the internal extraction control RNA directly to the unprocessed biological sample as this will lead to degradation and a loss in signal.

1. Add 4µl of the Internal extraction control RNA (BLUE) to each sample in RNA lysis/extraction buffer per sample.

2. Complete RNA extraction according to the manufacturer's protocols.

OneStep RT-qPCR detection protocol

For optimum performance and sensitivity.

All pipetting steps and experimental plate set up should be performed on ice. After the plate is poured proceed immediately to the OneStep amplification protocol. Prolonged incubation of reaction mixes at room temperature can lead to PCR artifacts that reduce the sensitivity of detection.

1. **For each RNA sample prepare a reaction mix according to the table below: Include sufficient reactions for positive and negative controls.**

Component	Volume
oasig OneStep or PrecisionPLUS OneStep 2X RT-qPCR Master Mix	10 µl
Mumps primer/probe mix (BROWN)	1 µl
Internal extraction control primer/probe mix (BROWN)	1 µl
RNase/DNase free water (WHITE)	3 µl
Final Volume	15 µl

2. **For each RNA sample prepare an endogenous control reaction according to the table below (optional):**
This control reaction will provide crucial information regarding the quality of the biological sample.

Component	Volume
oasig OneStep or PrecisionPLUS OneStep 2X RT-qPCR Master Mix	10 µl
Endogenous control primer/probe mix (BROWN)	1 µl
RNase/DNase free water (WHITE)	4 µl
Final Volume	15 µl

3. **Pipette 15µl of these mixes into each well according to your qPCR experimental plate set up.**
4. **Pipette 5µl of RNA template into each well, according to your experimental plate set up.**
For negative control wells use 5µl of RNase/DNase free water. The final volume in each well is 20µl.

5. If a standard curve is included for quantitative analysis prepare a reaction mix according to the table below:

Component	Volume
oasig OneStep or PrecisionPLUS OneStep 2X RT-qPCR Master Mix	10 µl
Mumps primer/probe mix (BROWN)	1 µl
RNase/DNase free water (WHITE)	4 µl
Final Volume	15 µl

6. Preparation of standard curve dilution series.

- 1) Pipette 90µl of template preparation buffer into 5 tubes and label 2-6
- 2) Pipette 10µl of Positive Control Template (RED) into tube 2
- 3) Vortex thoroughly
- 4) Change pipette tip and pipette 10 µl from tube 2 into tube 3
- 5) Vortex thoroughly

Repeat steps 4 and 5 to complete the dilution series

Standard Curve	Copy Number
Tube 1 Positive control (RED)	2 x 10 ⁵ per µl
Tube 2	2 x 10 ⁴ per µl
Tube 3	2 x 10 ³ per µl
Tube 4	2 x 10 ² per µl
Tube 5	20 per µl
Tube 6	2 per µl

7. Pipette 5µl of standard template into each well for the standard curve according to your plate set-up

The final volume in each well is 20µl.

OneStep RT-qPCR Amplification Protocol

Amplification conditions using oasig OneStep or PrecisionPLUS OneStep 2X RT-qPCR Master Mix.

	Step	Time	Temp
	Reverse Transcription	10 min	55 °C
	Enzyme activation	2 min	95 °C
Cycling x50	Denaturation	10 s	95 °C
	DATA COLLECTION *	60 s	60 °C


* Fluorogenic data should be collected during this step through the FAM and VIC channels

Interpretation of results


Target (FAM)	Internal control (VIC)	Positive control	Negative control	Interpretation
≤ 30	+ / -	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
> 30	+	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
> 30	-	+	-	POSITIVE QUALITATIVE RESULT do not report copy number as this may be due to poor sample extraction
-	+	+	-	NEGATIVE RESULT
+ / -	+ / -	+	≤ 35	EXPERIMENT FAILED due to test contamination
+ / -	+ / -	+	> 35	*
-	-	+	-	SAMPLE PREPARATION FAILED
+ / -	+ / -	-	+ / -	EXPERIMENT FAILED

Positive control template (**RED**) is expected to amplify between Cq 16 and 23. Failure to satisfy this quality control criterion is a strong indication that the experiment has been compromised.

*Where the test sample is positive and the negative control is positive with a Cq > 35, the sample must be reinterpreted based on the relative signal strength of the two results:

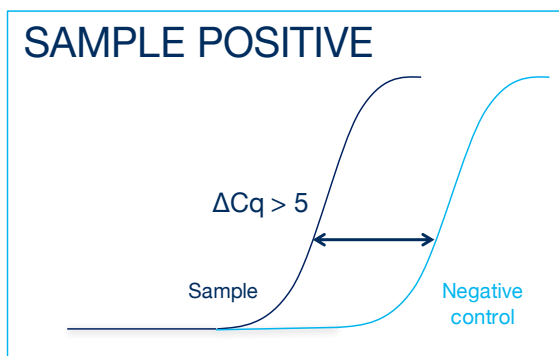

If the sample amplifies > 5 Cq earlier than the negative control then the sample should be reinterpreted (via the table above) with the negative control verified as negative.

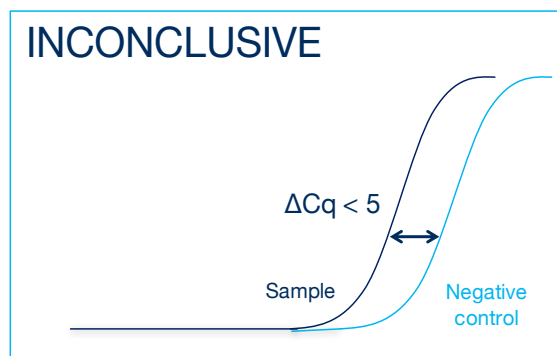

If the sample amplifies < 5 Cq earlier than the negative control then the positive sample result is invalidated and the result should be determined inconclusive due to test contamination. The test for this sample should be repeated.

Internal PCR control

The Cq value obtained with the internal control will vary significantly depending on the extraction efficiency, the quantity of RNA added to the RT and PCR reaction and the individual machine settings. Cq values of 28 ± 3 are within the normal range. When amplifying a Mumps sample with a high genome copy number, the internal extraction control may not produce an amplification plot. This does not invalidate the test and should be interpreted as a positive experimental result.

Endogenous control

The signal obtained from the endogenous control primer and probe set will vary according to the amount of biological material present in a given sample. An early signal indicates the presence of a good yield of biological material. A late signal suggests that little biological material is present in the sample.